
Instruction — LEAF & MAPS Adjustment

Freedom Designs, Incorporated
7801 S. Jackson Rd., Pharr, Texas 78577
phone: 800-331-8551 fax: 888-582-1509
e: customerservice@freedomdesigns.com
www.freedomdesigns.com

Tools needed:
- 5/32" Hex key wrench
- 3/16" Hex key wrench
- 7/16" Combination wrench
- 1/2" Combination wrench

LEAF Adjustment Instructions:

1. Loosen center (1) 1/4-20 socket head screw using a 3/16" Allen wrench (do not remove).
2. Remove rear 1/4-20 socket head screw using the 3/16" Allen wrench and re-position for left, right or midline position as shown above. Tighten both screws securely. (Note: for midline position, rear screw must be moved to the forward position.

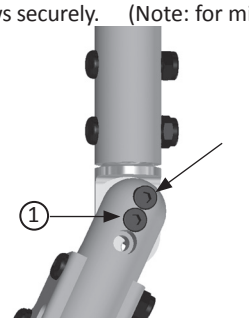

Rear screw position for frame set in abduction or windswept



Rear screw position for frame in midline





Abduction




Windswept Right




Windswept Left

MAPS Adjustment Instructions:

1. To adjust footplate angle position, loosen nuts (1) using a 1/2" combination wrench. Position footplate and re-tighten.
2. To adjust footplate depth, loosen nuts (2) using a 7/16" combination wrench while holding bolts (3) using a 5/32" Allen wrench.

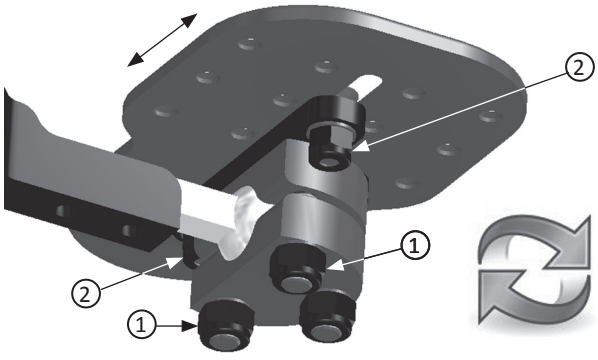

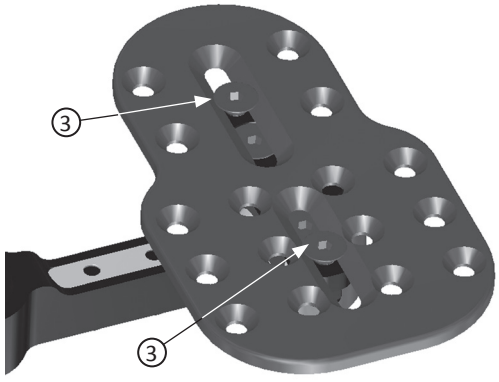



Eversion



Inversion



Rotation



Plantar / Dorsi Flexion